

SKIN RESURFACING

CO₂ Fractional laser

Laser Beauty Machine



FUNCTION



Fractional CO² Laser System

- ✓ Skin renewal and regeneration
- ✓ All colors tattoo removal
- ✓ Remove all kinds of wrinkles
- ✓ Remove acne and scars
- ✓ Remove pigmentation
- ✓ Treat telangiectasis
- ✓ Skin renewing
- ✓ Actinic keratoses
- ✓ Face lift, tighten and whiten skin
- ✓ Pierce ears
- ✓ Skin hydroponic, skin rough
- ✓ Teeth whitening
- ✓ Vaginal tightening
- ✓ E.N.T treatment




PARAMETER

Screen	10.4 inch LCD real color touch screen
Light conducting system	7 jointed arms to conduct light
Available treatment head	Scan、UltraPulse、Gynecology、Andrology
Wavelength	10.6 μ m
Power	40W
MODE	Scan, Normal
Laser energy	1–400MJ
Step value	1mj
Scan pattern size	0.1*0.1MM–20*20MM
Dot quantity	40000 dots at most
Pulse duration	0.2–4ms(Fractional);0.1~10ms(Normal)
Spot distance	0.1–2mm
Pulse energy	1mj to 400mj
Cooling system	air cooling system
Laser apparatus	Sealed off laser device stimulated by direct current
Condenser focus	f=100mm
Spot size	2.0 $\pm$ 0.2mm
Radiation time	0.011s
Main machine size	45*33*113cm
Net weight	45*33*113cm

ADVANTAGE

1 7 joint (Korea) make torsional spring light guide arm could ensure the effect of laser accurately treatment.

2 USA Coherent laser device, bigger and more stable power, longer life.

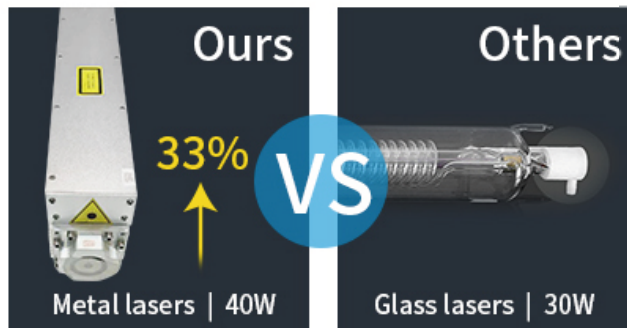

3 4 function: Lattice laser mode & pulse laser mode & Private mode & Vulvar treatment.

4 Air pump system: the air blowpipe puff away the smoke so operator's vision will more better.

5 The lattices laser mode we original use galvanometer. It can deal with: deep skin off & Tiny-exfoliative & non-exfoliative.



TECHNICAL FEATURES



Exclusive scanning heads for different skin problems. Minimum wounds. Perfect target the skin problem.

7-jointed arm light conduct, more convenient operating, less costs. 12 shapes for chosen. Support for hand-drawn free graphics.

Continuous, super pulse, fractional different Operation Modes, There is a wider use in clinical application.

Vaporization, ablation type, there is the gasification. Such as acne scar; stretch mark.

MICRO-EXFOLIATION LATTICE LASER:

The stratum corneum is intact, and the epidermis layer below it has a gasification. Such as: tattoo removal.

NON-EXFOLIATING LATTICE LASER:

Non-vaporized, non-ablative, with micro-thermal damage, no gasification. Such as: interface spot.



APPLICATION



6/9



COMPANY PROFILE



Ample spare parts reserve
Rich shell options
Instant delivery time



UNT has a strong research and development team, including electronics engineers (designer; AVR programming; circuit layout) construction engineer, UI designer engineer, QC control engineer...etc. We provide OEM/ODM design and program update.



In UNT we have been at the forefront of laser diode technology since 2010, We are committed with the clinics and practices to deliver them a platform that gives them the ability to get there results that the customers demand.

MORE PRODUCTS





Beijing UNT Technology CO.,Ltd

Address: No.M2-1, Xinggu Development Zone, Pinggu District, Beijing China. 101200

Website: <https://untlaser.en.alibaba.com/>

www.untlaser.com

Contact: 86-010-61960002

